


ENT

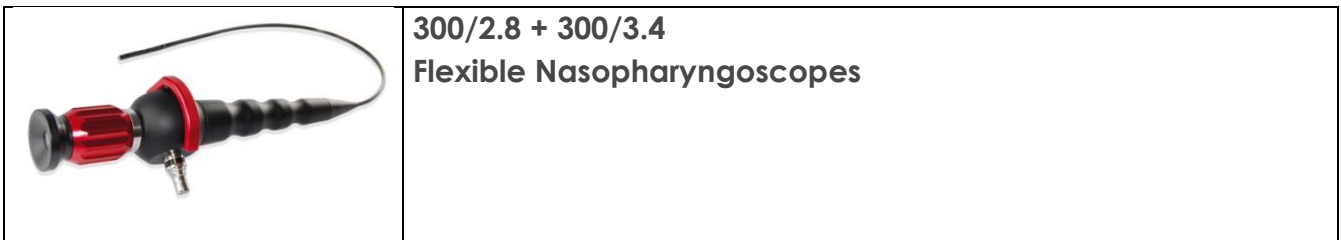


The fiber optic nasopharyngoscope offers optimal image quality by a multi-band optical system of highest precision (made in Germany). The extremely thin and flexible insertion tube ensures precise working in inaccessible regions. Wide angle of view with very good depth of field. Uniform illumination thanks to four light guide branches. Balanced control unit with ergonomic design. Standardized leak proof connector.

- >> **Absolutely leak proof**
- >> **Brilliant image**
- >> **Corrosion-resistant**
- >> **Ergonomic design**
- >> **18.000 pixel**
- >> **Available in the colors: red – gold – silver – blue**

Working length	300 mm	300 mm
Diameter	2,8 mm	3,4 mm
Angle of view	70°	70°
Bending	up/down 150°	up/down 150°
Article no:	50 28 30	50 34 30

Althea Deutschland GmbH

Im Goldäcker 14 | D-88630 Pfullendorf | Tel. +49 (0) 7552 92801-0 | Fax. +49 (0) 7552 92801-11
E-Mail: de.info@althea-group.com | www.althea-group.com/de

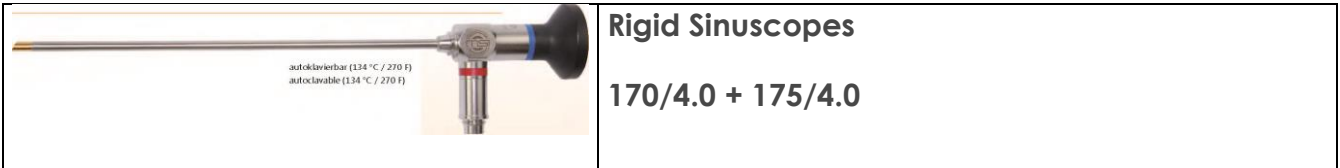

Rigid Sinuscopes
170/4.0 + 175/4.0

Completely made of high quality steel for highest corrosion resistance. All mechanical bonds are laser welded. Distal and proximal cover lenses consist of sapphire glass and are soldered in. That means: no adhesive surfaces, no leakage and no scratching. "High-resolution" optic systems for special brilliant image quality. Althea endoscopes are compatible with shafts of all different manufacturers and can be integrated in your present equipment / accessories without any problems. Different lines of vision are marked in terms of color. The blue ring shows that the endoscope is autoclavable.

- >> **Cover lens made of sapphire glass**
- >> **Corrosion-resistant**
- >> **134°C autoclavable**
- >> **Brilliant image**
- >> **Absolutely leak proof**

Compatible	STORZ	STORZ	STORZ	STORZ	STORZ	STORZ	STORZ
Diameter	4,00 mm	4,00 mm	4,00 mm	4,00 mm	2,70 mm	2,70 mm	2,70 mm
Working length	175 mm	175 mm	175 mm	175 mm	187 mm	187 mm	175 mm
Angle of view	0°	30°	45°	70°	0°	30°	0°
Aperture angle	95°	95°	95°	95°	80°	80°	80°
Article no:	10 40 00	10 40 30	10 40 45	10 40 70	10 28 00	10 28 30	11 28 00

Compatible	STORZ	STORZ	STORZ	STORZ
Diameter	2,70 mm	2,70 mm	2,70 mm	2,70 mm
Working length	175 mm	110 mm	110 mm	74 mm
Angle of view	30°	0°	30°	0°
Aperture angle	80°	80°	80°	80°
Article no:	10 28 30	12 28 00	12 28 30	50 28 00

Althea Deutschland GmbH

Im Goldäcker 14 | D-88630 Pfullendorf | Tel. +49 (0) 7552 92801-0 | Fax. +49 (0) 7552 92801-11
E-Mail: de.info@althea-group.com | www.althea-group.com/de